



Standard Colors:



Velvet Red
(Installed)



Coastal Blue
(Included in Box)

NEW Sporty Optional Colors:




Crisp White



Canary Yellow




Stone Grey

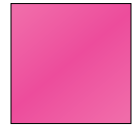

Orchid Pink



Army Green
(Matte)

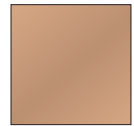

Desert Tan

Buzzaround GB120



The Buzzaround CarryOn fold-flat scooter is small in stature and big on fun! This super portable, sporty scooter will take you where you want to go in comfort and style! It features a 300 lb. weight capacity, an infinite adjustable tiller with a convenient USB port and drink holder, full front and rear light package, and flip up armrests. Enjoy long range coverage up to 14 miles on the airline friendly 12.5 AH lithium-ion battery, which can be removed and taken as a carry on while traveling by air, sea or train in a stylish carry case. A new rear basket is mounted on the deck to hold an optional spare 12.5 AH battery for even longer range and serves as a kickstand for storing upright. The heaviest piece is 47.5 lbs. (without seat and batteries) and the height when folded is only 12". It disassembles easily into two pieces and is designed to store in the smallest compartments. It comes standard with red color shrouds installed and an extra set in blue. New optional color shroud kits let you express a style all your own! Pack your dreams and go anywhere with the CarryOn Scooter!

**NEW & IMPROVED
Sporty Fold-Flat Scooter!
Airline Friendly!
USB Port Standard!
Optional Travel Bag!**

Standard Features:

- 4-point carry handles for easy transport
- 12" height when folded (without seat)
- 4" ground clearance
- 300 lb. weight capacity
- Flip up armrests
- Infinite adjustable tiller
- USB port & drink holder
- Disassembles easily
- Up to 14 mile range
- Easy access, airline friendly 12.5 AH lithium ion battery pack





LIGHTWEIGHT PORTABLE SCOOTER

Model	GB120
Medicare Code	Not Coded
Weight Capacity	300 lbs
Drive Wheels	Rear
Maximum Speed	3.6 mph
Operating Range ¹	14 miles
Ground Clearance Under Deck	4"
Ground Clearance to Transaxle/Motor	1.5"
Turning Radius	48.75"
Battery Type (Lithium Ion)	12.5 AH
Colors	Red, Blue, White, Yellow, Grey, Pink, Green & Tan
Freewheel Mode	Yes
Electronic Speed Control	Yes
Electro-Mechanical Brakes	Yes
Charger	Off Board
Controller	50 Amp
Dimensions	
Length	38"
Width	19.5"
Height (ground to top of back of seat)	34.5"
Maximum Height (folded without seat)	12"
Ground to Top of Seat	23.5"
Ground to Top of Deck	6"
Deck to Top of Seat	17.5"
Front Axle to Rear Axle	28.7"
Weight of Unit: (assembled)	65.5 lbs
Main Body	47.5 lbs
Seat w/Arms	12 lbs
Batteries	12.5 AH 6 lbs
Tires:	Solid/Low Profile
Front	6"
Rear	7"
Rear Anti-Tip Wheels	2"
Standard Seat	Stadium Style
Back Height - no headrest	13.25"
Width x Depth	15" x 13.75"
Color	Black Vinyl
Warranty:	
Frame	Lifetime
Drive Train	2 Years
Electronics ²	2 Years
Batteries	1 Year

¹ Battery range at 200 lbs., and will vary due to rider weight, drive surface, terrain, and battery type.
² Electronics warranty excludes batteries.



Convenient kickstand for small space storage



Folds easily for a compact profile of only 12" high (without seat)



Optional spare battery storage in new rear basket



Front & rear lighting package



Front carry handle makes travel easy



Conveniently located USB port



LED battery indicator on the console



NEW! Optional Travel Bag



Order the optional Golden Travel Case for on the go storage!

BUZZ GB120



Scan for Video

Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.



Phone: 800-624-6374
 Fax: 800-628-5165
 www.goldentech.com

Find us on Facebook
www.facebook.com/goldentechnologies

SHT-GB120-081324