



Standard Colors



Grey



Teal

**Remote
Control App
for Enhanced
Convenience**

**Folds Easily
Drives Remotely
Tight Turning
Radius**

Standard Features:

- Smart Control Mobile App
- Memory Foam Seat Cushion
- Padded Armrests
- Adjustable Recline Back
- Spacious Leg & Foot Area
- Length Adjustable Footrest
- Front Suspension
- Omni-Directional Wheels
- Airline Friendly, Lithium-Ion Battery
- Easy-to-Reach Storage Compartment
- System Monitoring Voice Alerts
- USB Charging Ports
- Front & Rear Lighting
- Water Resistant Backpack

**Ally Pro
GP304**



Visit App Store
for Remote
Control



The Golden Ally Pro™ is next in line for our smart power wheelchair series. It features a modern design that applies robotics technology for easy folding and remote control. Upgrades from the Ally™ include 10" omni-directional front wheels allowing for more leg and foot room, the ability to climb larger obstacles, and rolling forward and sideways enabling the chair to turn in a tight radius. It offers great ride quality and capability, compliments of its large, pneumatic rear drive wheels and shock-absorbing front suspension. The two-position adjustable back angle provides long-lasting comfort. Arms are adjustable for additional seat width and feature padded armrests. Its comfortable seating and foot positions will also be a welcome change from the compromises of the typical power chair. The Ally Pro is perfect for small spaces indoors or open spaces outside. With a long-lasting lithium-ion battery, user-friendly interface, and Smart Control Mobile App, you can navigate any environment with ease and confidence. Move effortlessly with the Ally Pro, your perfect travel partner!

**Exclusive
Omni-Directional
Wheels with
Self-Cleaning
V-Groove!**



Model Number	GP304
Medicare Code	Not Coded
Weight Capacity	330 lbs.
Drive Wheels	Rear
Maximum Speed	5 mph
Operating Range ¹	13 miles
Ground Clearance	2.5"
Turning Radius	37"
Type Batteries	12 AH Lithium-Ion
Colors	Teal, Grey
Freewheel Mode	Yes
Electronic Speed Control	Yes
Electro-Mechanical Brakes	Yes
Charger	Off Board
Controller	40 Amp Remote
Dimensions:	
Length	39"
Width	25"
Width Between Arms	20" max.
Height (ground to top of back of seat)	37"
Ground to Top of Seat	20.5"
Weight of Unit: (assembled)	75 lbs.
Heaviest Piece	69.6 lbs.
Battery	5 lbs.
Tires:	
Drive Wheels	10" Pneumatic
Casters	10" Omni-Directional
Anti-Tip Wheels	2"
Standard Seat:	
Back Height - no headrest	19"
Width x Depth	17" x 16"
Color	Black Fabric
Footrest:	
Height adjustable	Yes
Warranty:	
Frame	5 Years
Drive Train	13 Months
Electronics ²	13 Months
Batteries	12 Months

¹ Battery range tested at 200 lbs. and will vary due to rider weight, drive surface, terrain and battery type.


² Electronics warranty excludes batteries.

FDA Class II Medical Devices designed by Golden Technologies are intended to assist persons with physical disabilities.

Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.



Padded Armrests

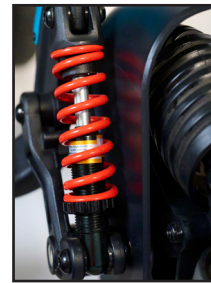


Front Suspension



Highly Visible Rear Lighting



Zippered Water Resistant Backpack

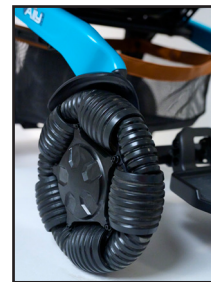

Omni-Directional Wheels with Self-Cleaning V-Groove



Pre-wired Controller for right or left side placement

Optional Accessories:



Extendable Headrest



Cup/Phone Holder



Heated Seat Control



Wireless Remote



Phone: 800-624-6374

Fax: 800-628-5165

www.goldentech.com



www.facebook.com/goldentechnologies

LIT-SHT-GP304, Rev.A 012425